


Dependa-Bar™

A revolutionary product which combines the sturdiness of a wall mounted grab bar and the functionality of a dual support rail that moves in step during standing and transferring in the bathroom.

Features and Benefits

- **Ease of Installation** – Mounts using conventional means to walls of most constructions. Bolt pattern compatible with Winglts® Fastening System (sold separately) if wall stud location is not appropriate.
- **Pivoting & Locking Technology** – Dependa-Bar™ locks every 45° (5 positions), providing a support rail that moves with user. When not in use, the Dependa-Bar™ conveniently locks against the wall. Ideal where bathtub is used by family members who do not need the support.
- **Versatile** – Can be installed in the bathtub or beside the toilet.
- **Water Resistant** – Stainless steel construction provides an easy to clean, durable grab rail that provides years of life.
- **Safe** – Fixed horizontal grab bar has ADA compliant diameter (1½”) and wall clearance (1½”). When locked against the wall, the Dependa-Bar™ is aligned with the fixed horizontal grab bar.

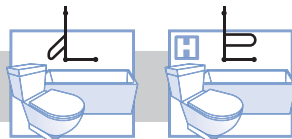

Bathroom

Hospital



Technical Specifications:

Pivot Rail	1.25" / 3.175 cm grip diameter, #304 Stainless steel construction.																		
Pivot Rail Height	30½" / 77.74 cm vertical rail height.																		
Pivot Rail Length	"P" shaped dual rail extends 18" / 45.72 cm length from bathtub wall (6" / 15.24 cm spacing between horizontal rails). Locks every 45° (5 positions).																		
Horizontal Grab Bar	1.25" / 3.175 cm grip diameter. 1.5" / 3.81 cm clearance from wall (ADA compliant). Height adjustable at time of installation, typically 4" / 10.16 cm from top of bathtub surface. #304 Stainless steel construction.																		
Grab Bar Lengths	18" / 45.72 cm and 32" / 81.28 cm lengths available for lower grab rail.																		
Weight Capacity	250lbs / 113 kg.																		
Assistance Level	<table border="1"> <thead> <tr> <th>M</th> <th colspan="4">ASSISTANCE LEVEL</th> <th>M</th> </tr> </thead> <tbody> <tr> <td>I</td> <td>1</td> <td>2</td> <td>3</td> <td>4</td> <td>A</td> </tr> <tr> <td>N</td> <td></td> <td></td> <td></td> <td></td> <td>X</td> </tr> </tbody> </table>	M	ASSISTANCE LEVEL				M	I	1	2	3	4	A	N					X
M	ASSISTANCE LEVEL				M														
I	1	2	3	4	A														
N					X														

*Patents Pending.