


# SuperPole™ with SuperBar™

A horizontal pivoting and locking bar that mounts to the SuperPole™, suitable for users who require moderate assistance to achieve small steps during transfers.

## Features and Benefits

- **Support That Moves With You** - To allow for 2 and 3 step transfers, the SuperBar™ pivots and locks into 8 different positions around the SuperPole™. It provides excellent stability while in the locked position. To allow for support throughout the transfer process, pull laterally to assist with standing and sitting, or lift, pivot and lock into the next position. The SuperBar™ automatically locks when released to prevent forward or lateral falls.
- **Support Where It's Needed** - The unique SuperPole™ with SuperBar™ combination provides gripping, leaning and weight bearing support in front of the user, exactly where it's needed. Well located support helps to reduce the risk of fall related injuries by eliminating the need to reach, twist or sometimes lunge for support.
- **Reduce Caregiver Strain** - Allows the caregiver to supervise a transfer with minimal physical exertion.

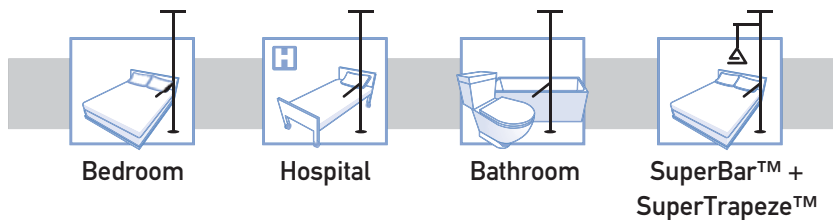

## Technical Specifications:

SuperBar™	1.5" / 3.81 cm grip diameter, 16" / 40 cm length, locks every 45° (8 positions), Height adjustable, typically 24"-38" / 61-97 cm high with 96" ceilings.
Floor to Ceiling Range	93" to 99" / 236 to 253 cm (standard).
Uni-Fit Extender™	Increase floor to ceiling range from 100" to 120" / 254 to 305 cm.
SuperPole™ Ultra	NOT TO BE USED WITH THE SUPERBAR™.
Ceiling Plate Extender™	Allows standard ceiling plate 20" / 51 cm length to span 30" / 76 cm with holes for ceiling attachment at 24" / 61 cm ceiling beam centers.
Pole Material	Heavy gauge steel, 1.5" / 3.81 cm diameter. High temperature baked powder coat gloss white enamel.
Grip Material	Grey, premium quality, 20" / 51 cm long. Does not contain natural dried rubber or latex.
Expansion Mechanism	Steel jackscrew, zinc plated, 1-1/8" fine cut thread.
Weight Capacity	250lbs / 113 kg. Heavy Duty weight capacity 450 lbs / 205 kg.
Assistance Level	