



SuperPole™ with SuperTrapeze™

An offset trapeze handle that mounts to the SuperPole™. Suitable for users who require moderate assistance to achieve a sitting position in bed, with all the benefits of the SuperPole™.

MEDICARE BILLING
CODE #E0940
(Freestanding trapeze)

Features and Benefits

- **Innovative** - The “e2”™ trapeze handle provides a second hand rail above a base rail, which simulates two rungs in a ladder. This provides better leverage for those individuals with limited arm strength to more easily “climb up” to a seated position, at the same time reducing the energy and range of motion normally required when using a standard trapeze bar.
- **Storage** - The open design of the trapeze handle allows it to be hooked out of the way onto the SuperPole™ when not in use.
- **Space Saver** - Unlike conventional trapeze systems, the SuperPole™ with SuperTrapeze™ installs beside the bed and does not force the bed out from the wall, saving valuable space in smaller rooms.
- The “e2” trapeze handle is also sold separately with an adjustable strap and hook for upgrading or replacing existing trapeze handles.

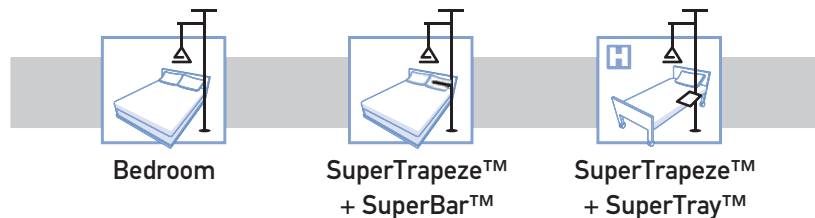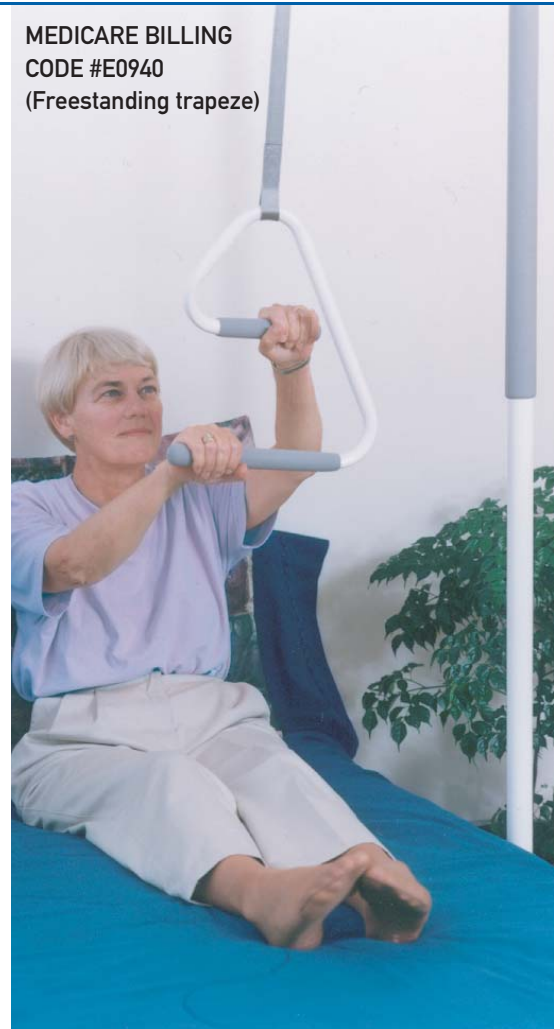

Technical Specifications:

| | | | | | | | | | | | | | | | | | | | |
|-------------------------|--|---|------------------|---|---|--|---|---|---|---|---|---|---|---|--|--|--|--|---|
| SuperTrapeze™ | 0.875" / 2.2 cm diameter, Adjustable 1" / 2.5 cm nylon strap attaches to welded support arm (adjustable in height on SuperPole™). | | | | | | | | | | | | | | | | | | |
| Floor to Ceiling Range | 93" to 99" / 236 to 253 cm (standard). | | | | | | | | | | | | | | | | | | |
| Uni-Fit Extender™ | NOT TO BE USED WITH SUPERTRAPEZE™. | | | | | | | | | | | | | | | | | | |
| SuperPole™ Ultra | NOT TO BE USED WITH SUPERTRAPEZE™. | | | | | | | | | | | | | | | | | | |
| Ceiling Plate Extender™ | Allows standard ceiling plate 20" / 51 cm length to span 30" / 76 cm with holes for ceiling attachment at 24" / 61 cm ceiling beam centers. | | | | | | | | | | | | | | | | | | |
| Pole Material | Heavy gauge steel, 1.5" / 3.81 cm diameter. High temperature baked powder coat gloss white enamel. | | | | | | | | | | | | | | | | | | |
| Grip Material | Grey, premium quality, 20" / 51 cm long. Does not contain natural dried rubber or latex. | | | | | | | | | | | | | | | | | | |
| Expansion Mechanism | Steel jackscrew, zinc plated, 1-1/8" fine cut thread. | | | | | | | | | | | | | | | | | | |
| Weight Capacity | 250lbs / 113 kg. Applied weight of 150 lbs / 68 kg. | | | | | | | | | | | | | | | | | | |
| Assistance Level | <table border="1"> <tr> <td>M</td> <td colspan="4">ASSISTANCE LEVEL</td> <td>M</td> </tr> <tr> <td>I</td> <td>1</td> <td>2</td> <td>3</td> <td>4</td> <td>A</td> </tr> <tr> <td>N</td> <td></td> <td></td> <td></td> <td></td> <td>X</td> </tr> </table> | M | ASSISTANCE LEVEL | | | | M | I | 1 | 2 | 3 | 4 | A | N | | | | | X |
| M | ASSISTANCE LEVEL | | | | M | | | | | | | | | | | | | | |
| I | 1 | 2 | 3 | 4 | A | | | | | | | | | | | | | | |
| N | | | | | X | | | | | | | | | | | | | | |